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**Instructions for use and handling**

1. Use the Derma Roller as described and avoid using it in the following situations:

- Do not use on wounds, rashes, inflammation or other skin irritations.

- Keloid skin.

2. Wash your hands and the area to be treated thoroughly.

3. Apply only with light pressure.

4. Apply the roller only for a short time in the same areas.

5. The skin around lips and eyes is thinner, avoid using in these areas.

6. Do not use on children and animals.

7. If there are signs of rash, inflammation, allergies, etc., discontinue use and consult a doctor.

8. If cosmetic liquid gets into your eyes, wash them out thoroughly immediately.

9. Do not expose skin to excessive UV light after application and apply sunscreen.

10. Avoid sauna and other similar places after application.

**Storage and Handling Precautions**

1. Keep the cap closed and store the product in a dark place.

2. Keep the product out of the reach of children and animals.

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**Use:**

Place the roller attachment on the active substance container. Repeated pressing on the attachment pumps active substance onto the roller. Then the active substance can be applied with the roller to stimulate the skin.

**Roller attachment information**

Number of needles:

540 pins

Length of needles:

0.25 mm

Purpose of use: Absorption

cosmetics and skin

stimulation.

**Serum information**

Content: 10 ml / 0.34 fl

oz Use by: Within 2 years

from the date of production

and within 1 month

after opening Functional

Cosmetics-Brightening

Functional Cosmetics, Anti-

wrinkle Extract Serum

Effective: brightening and

wrinkle improvement.

**How to use:**

Apply the cosmetics to the face to accelerate absorption.

**Ingredients:**

Water, Niacinamide, Hydroxyethylcellulose, Chlorphenesin, Usnea Barbata (Lichen) Extract, Zanthoxylum Piperitum Fruit Extract, Pulsatilla Koreana Extract, 1,2-Hexanediol, Madecassic Acid, Asiatic Acid, Hydrogenated Phosphatidylcholine, Sucrose stearate, Sodium Hyaluronate, DNA, Butylene Glycol, Glycerin, Panax Ginseng Root Extract, Beta-Glucan, Centella Asiatica Extract, Pentylene Glycol, Polyglutamic Acid, Serine, Madecassoside, sr-Spider Polypeptide-1, Propolis Extract, Aspartic Acid, Acetyl Hexapeptide-8, Adenosine, Leucine, Alanine, Lysine, Arginine, Tyrosine, Phenylalanine, Valine, Threonine, Proline, Isoleucine, Histidine, Methionine, Cysteine, Polysorbate 20, Lecithin

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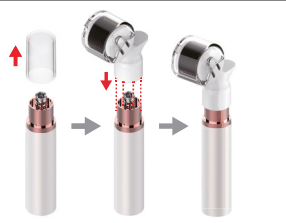
Unpack the roller attachment and active substance container.

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Remove the cap and place the roller attachment on the active substance container.

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To avoid infection, disinfect or sterilize the roller attachment. A cotton pad soaked in ethanol or alcohol is best suited.

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Separate the roller attachment from the active substance container.

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Wash the inside of the roller attachment (connecting piece to the active substance container) with warm water.

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After drying, close the protective cap of the roller and the plastic cap of the serum and store them in a cool place out of direct sunlight.

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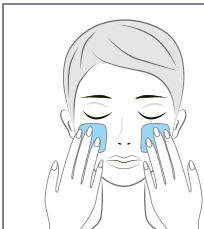
Wash the face and then apply a toner in preparation.

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
Remove the protective cap from the roller attachment and disinfect it with alcohol or ethanol.

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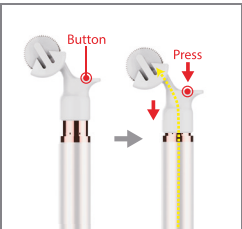
Place the roller attachment on the active substance container and pump 3-5 times so that the active substance flows onto the roller.

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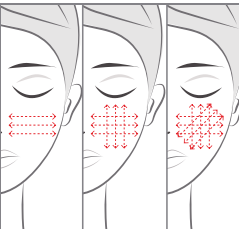
Apply the active ingredient in horizontal, vertical and diagonal movements.

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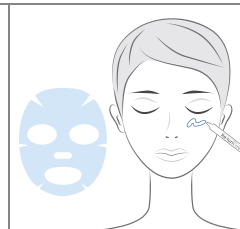
After the treatment, apply a collagen-rich cream or mask to soothe the skin.

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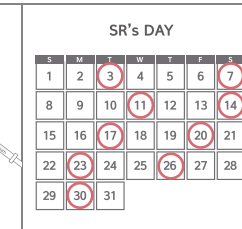Apply the Derma Roller 1-2 times a week for 10 minutes. Give the skin a break after 1 month of use.