
**TITANIA®**

REG. GERMAN BRAND

# THE BEAUTY FACER®


**BEDIENUNGSANLEITUNG**  
**OPERATING INSTRUCTIONS**  
**MANUEL D'UTILISATION**  
**INSTRUCCIONES DE USO**  
**UPUTSTVO ZA UPOTREBU**  
**РУКОВОДСТВО ПО ЭКСПЛУАТАЦИИ**  
**دليل الاستخدام**

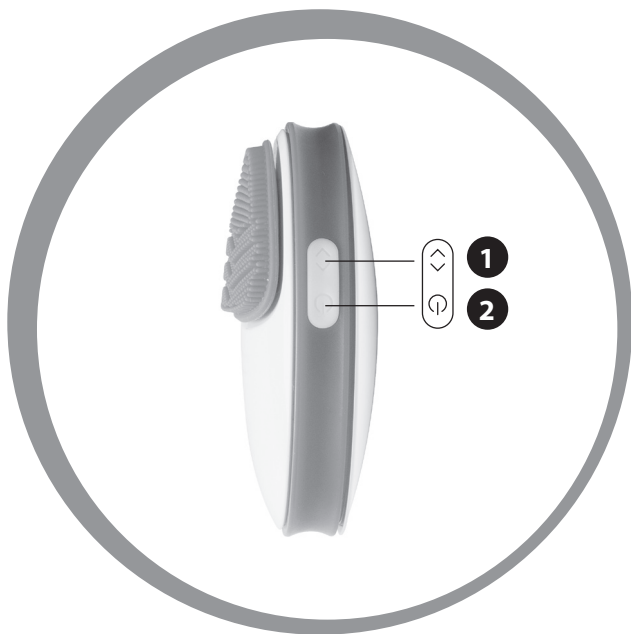- 1** 2 Geschwindigkeiten  
2 speeds  
2 vitesses  
2 velocidades  
2 brzine  
2 скорости  
سرعتان

- 2** An- / Aus-Schalter  
On/Off switch  
Interrupteur Marche/Arrêt  
Interruptor de encendido / apagado  
Prekidač za uljučivanje i isključivanje  
Переключатель вкл./выкл.  
مفتاح التشغيل/الإيقاف

Tägliche Reinigung und Anti-Aging  
Daily cleansing and anti-ageing  
Nettoyage quotidien et anti-âge  
Limpieza y tratamiento anti-  
envejecimiento diarios  
Dnevno čišćenje i anti-aging  
Ежедневная очистка и  
антивозрастной уход  
التنظيف اليومي للوجه ومقاومة الشيخوخة

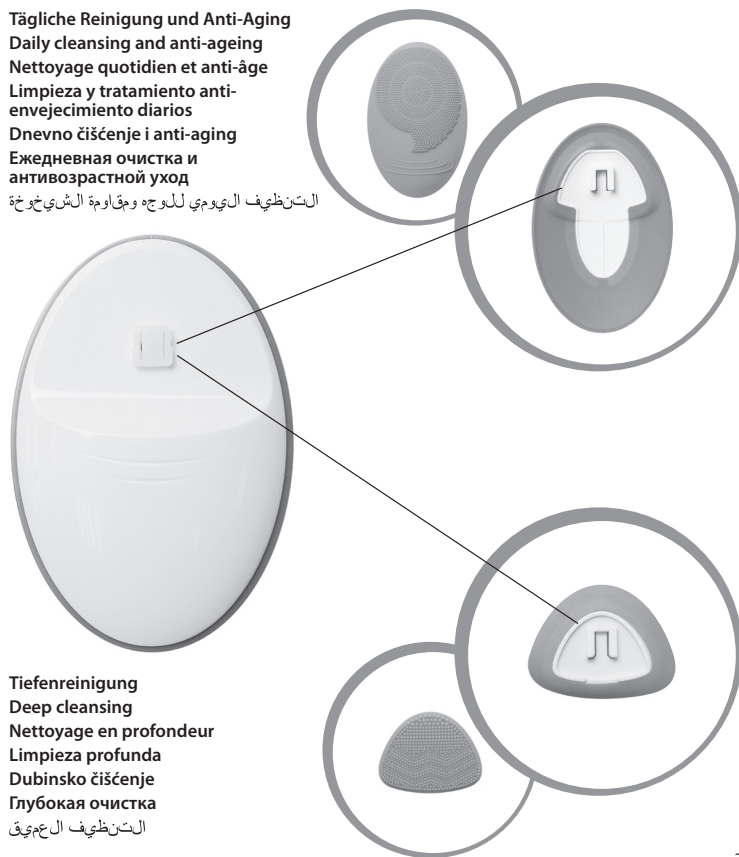Tiefenreinigung  
Deep cleansing  
Nettoyage en profondeur  
Limpieza profunda  
Dubinsko čišćenje  
Глубокая очистка  
التنظيف العميق

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6

## Warning (important notes)

Please read these operating instructions carefully and completely before using the product. Only use the product for the intended purpose as described in these operating instructions. These operating instructions must be stored and made easily available. The Beauty Facer® is not suitable for commercial, professional or medical applications, instead it is only designed for home use. The Beauty Facer® is not designed to be used outdoors.


## Careful/Caution!

Never use the device for more than 3 minutes. Avoid using cleansers with a silicone basis or granular cleansers since these may damage the silicone attachments. Avoid direct sunlight. Never expose the device to extreme heat or boiling water. Do not submerge in water. Keep away from children. Do not use on eyelids, avoid contact with eyes. Do not use on irritated or damaged skin. In the event of irritation, redness, itching or discomfort, switch the device off. In the event of skin disease, consult your doctor before using.

## Inserting and replacing batteries

Please check the type of battery to be inserted prior to use. Then please insert the batteries in accordance with the instructions depicted below.


1. Lift the rear flap to open the battery compartment and remove.
2. Insert the batteries into the compartment. Pay attention to ensure the batteries are the correct way round.
3. Replace the flap in such a way that the characters ||> are on top of one another and engage by pressing.

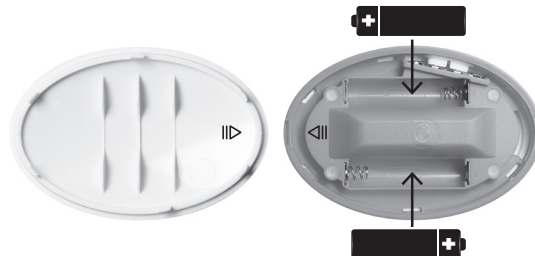7

## Use

1. Wash the face with warm water.
2. Apply your preferred facial care product to the silicone attachment.
3. Switch the device on using the **2** button .
4. Use the **1** button, to select one of the two speeds.
5. Position the device at the chin and start by making small circular motions upwards. Work along the jawline to the ear and then to the cheek. Repeat on the other side of the face and work across the forehead and along the nose.

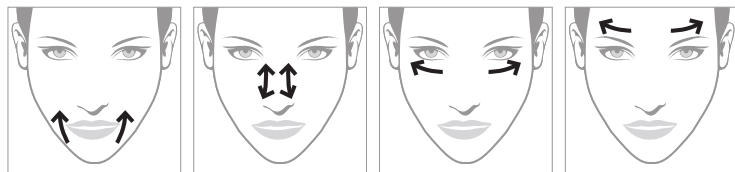6. Switch the device off using the **2** button.

For anti-ageing treatments, use the large silicone attachment on dry skin.

## Changing the silicone attachment

1. Using both thumbs, push the silicone attachment upwards, this will raise it out of the plastic spindle.
2. Place the other silicone attachment on the same spindle and push down with both thumbs.
3. To make cleaning easier, the silicone attachment can be removed from the plastic frame and then replaced prior to the next use.
4. Slide the large silicone attachment **[1]** onto the device for general, daily cleansing or a anti-ageing treatment.
5. Slide the small silicone attachment **[2]** onto the device for deep cleansing of the face.

## Cleaning


Do not use aggressive cleaning agents. Cleaning agents and soaps must be free of alcohol, benzene and acetone. Always switch the device off before cleaning. Clean the device after each use with dry or lightly moistened cloth, as required using water and a mild soap, then rinse. Make sure that no moisture is able to penetrate the device. Dry the device with a lint-free cloth/towel. The silicone attachments can be cleaned with a mild soap under running water once they have been removed from the device.

## Technical specifications

Battery-operated: 3 V, 2 x 1.5 V batteries of type AAA (included in the scope of delivery)

## Environmental Protection

Used devices with the depicted symbol may not be disposed of with the household waste. Used batteries and rechargeable batteries labelled with one of the depicted symbols may not be disposed of with the household waste. They must be disposed of via a collection point for used devices, used batteries or in the specialist waste (find out more from your community) or handed in at the retailer you originally purchased from. These places will ensure environmentally friendly disposal.
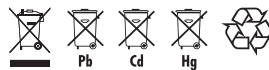
## Guarantee notes

The device has been subject to a careful final check. However, should you have a complaint despite this, please send the device back to us along with the purchase receipt. We offer a guarantee of 24 months as of the date of purchase. We do not accept liability for damages caused by incorrect handling, unintended use or wear. We reserve the right to make technical changes.

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**TITANIA\*** Fabrik GmbH